




- Etched resin-bonded bridge
- Minimal preparation
- Anterior lingual splints
- High strength of base alloy
- Simple pontic and two wing preparation

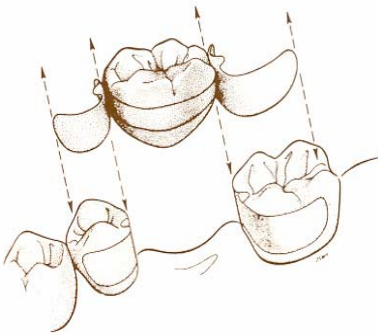

Prime consideration for design is to provide a distinct path of insertion. Design consists of the proximal and lingual segments and the occlusal rest.

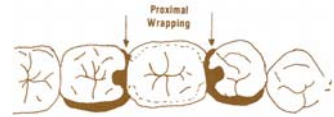

- 0.5 mm Lingual Reduction
- Proximal Notches

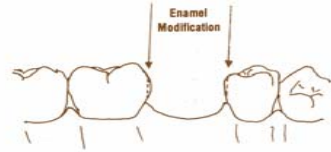

THE PROXIMAL SEGMENT:

Make a slight wrap of the proximal metal beyond the contact area toward the proximo-facial line angle. Develop a distinct vertical path of insertion.

THE LINGUAL SEGMENT:

The gingival edge of the lingual segment should be positioned supra-gingivally.

THE OCCLUSAL REST:

The posterior rests are somewhat smaller than those required for a removable partial denture. The main function is to dissipate the occlusal stresses across adjacent teeth.

Bonding:

Bond with Panavia Opaque (Panavia OP)

If you have any questions, call us

888-875-6011