


The preparation for a Procera® crown should follow the basic principles that apply for optimizing the strength and esthetics of an all-ceramic crown. **Chamfer prep is the optimum choice.**

- Anterior reduction should be 1mm—1.5 mm with rounded inner angles and edges.
- Posterior occlusal reduction should be 1.5 mm—2 mm and have a relatively flat surface topography.
- Procera is cemented with your material of choice.

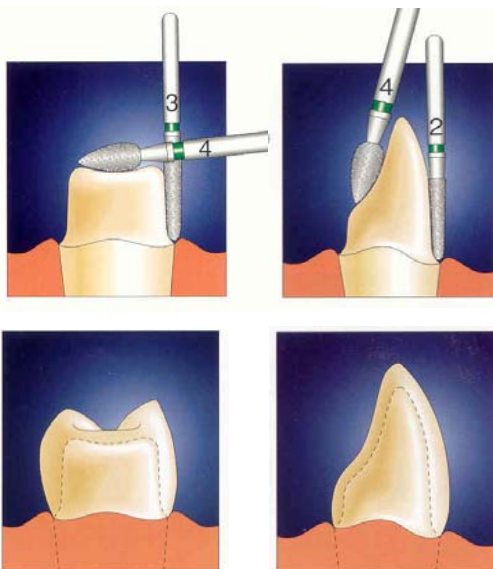

### **New Ultra Thin White, Ultra Translucent & Zirconia Copings**

New additions to the Regular Procera family are the ultra-thin 0.4 mm white, ultra-translucent 0.04 mm copings as well as Procera® Zirconia! Regular Procera copings are 0.6 mm.

Use Procera White coping Crowns for :

- Bleach Shades
- Cases with dark preps

Use Procera Translucent Crowns for:

- Teeth that have natural vibrant dentin
- Case that need more room for porcelain

Use Procera Zirconia Crowns for:

- Superior strength in the posterior

Use Regular Crowns for:

- All other applications

If you have any questions, call us

**888-875-6011**